

FRONT SURFACE

Cross-sectional view


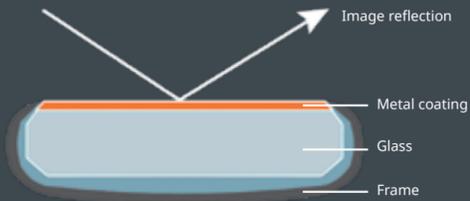


Image reflected directly on the surface of the glass, distortion-free
NO DOUBLE IMAGE

BACK SURFACE

Cross-sectional view of an aluminium mirror

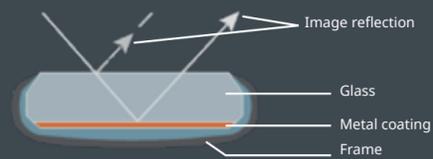

Reflecting surface under the glass
DOUBLE IMAGE



The finitions of the frame also matters

Flat and ergonomic soldering

More comfortable for the patient while retracting the cheek



Satined finition

Prevent light reflection in the eyes
Avoid eye fatigue



Optimal crimping

The frame suits the mirror
No sharp edges
No blood and dirt retention



www.acteongroup.com

Pictures:
1- Courtesy of Dr Tassery
2- Courtesy of Dr Cochet

PRODONT-HOLLIGER | A Company of ACTEON Group
La Marnasse | 63880 OLLIERGUES | FRANCE
Tel + 33 (0) 473 955 642 | Fax + 33 (0) 473 955 699
E-mail : info@prodont-holliger.fr | www.acteongroup.com

ACTEON

PUREREFLECT
SEE THE TRUE IMAGE

I AM
YOUR THIRD EYE

ACTEON

Non contractual document - Ref. D00125 - V1 - 08/2018 - Copyright © 2018 ACTEON®. All rights reserved. No information or part of this document may be reproduced or transmitted in any form without the prior permission of ACTEON®.

PURE REFLECT®
TECHNOLOGY
an exclusivity of ACTEON PRODONT-HOLLIGER

Inspired by the
aerospace industry

No double image (front surface)




Very luminous




No color and shape distortion
Accurate and sharp image

The authentic color reflection

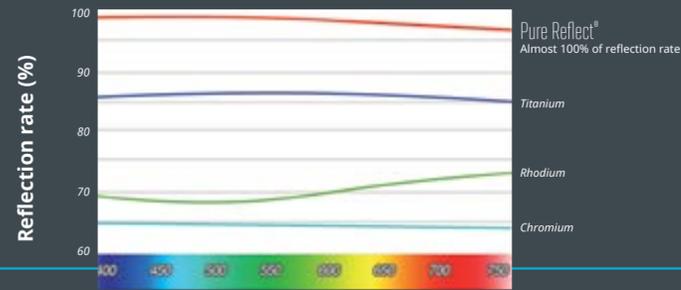

Visible wavelengths (nm)
*source: CNRS laboratory - France 12/2013

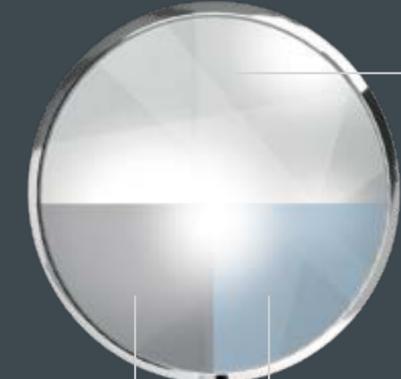


RHODIUM
a greyish reflection

TITANIUM
a blueish reflection

MOUTH MIRRORS

From diagnostic to treatment,
the benefit to all dental clinical acts

Visual examination is the first step to diagnostics
Deserve the true color reflection for an
optimal performance

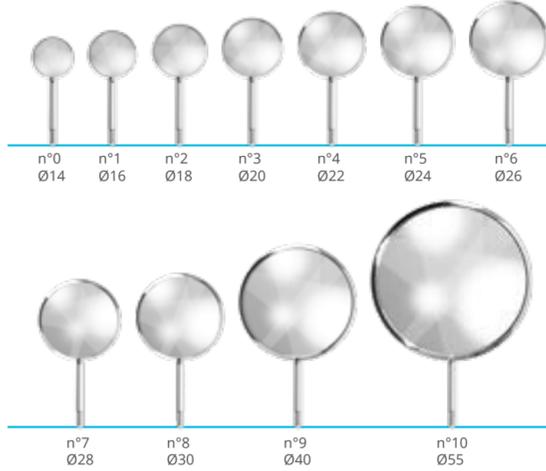

PURE REFLECT®



Competitor



1

MICROSCOPIC MIRRORS

The microscopic perfectly visible

- Make the access to the apex easier
- Appreciate the accurate vision of the pulp chamber
- Save space in the oral cavity: handle instruments and micro mirrors comfortably



Visibilité optimale
de l'isthme



2



Prof. H. Tassery, Marseille

"In my everyday clinical practice, the Pure Reflect® mirror remains a valuable tool permitting the perfect vision of the shades and details of dental structures. It permits working without any fatigue, with or without visual aids such as magnifying glasses and microscopes. The Pure Reflect® mirror has become essential and has made traditional mirrors lapsed."

PHOTOGRAPHIC MIRRORS




Prof. J.F. Lasserre, Bordeaux

"The exceptional luminosity and the colored fidelity define the Pure Reflect® mirrors as indispensable tools to a qualitative intrabuccal photograph"

Effective communication with patients

Make the difference
with the Pure
Reflect® technology

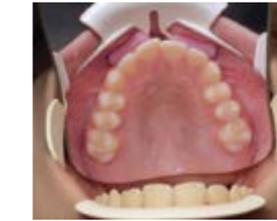

PURE REFLECT®



Titanium

- Highlight the dental treatment and restoration to promote your practice
- Take pictures of clinical cases for books, publications and patient records



Occlusal



Vestibular



Lingual



All our mirrors are double sided
Made out of glass
8 sizes available



Handle
Hold the photo mirror and
retract the cheek at the
same time